

Human sIL-4R ELISA KIT

Catalog No. 850.670.192

INTENDED USE

The Human sIL-4R ELISA is to be used for the in-vitro quantitative determination of soluble IL-4 receptor molecule (sCD124) in human serum, plasma, buffered solutions or cell culture medium. **This kit has been configured for research use only and is not to be used in diagnostic procedures.**

PRINCIPLE OF THE METHOD

The sCD124 Kit is a solid phase sandwich Enzyme Linked-Immuno-Sorbent Assay (ELISA). A monoclonal antibody specific for sCD124 has been coated onto the wells of the microtiter strips provided. Samples, including standards of known sCD124 concentrations, specimens and unknowns are pipetted into these wells.

During the first incubation, the sCD124 antigen and a biotinylated monoclonal antibody specific for CD124 are simultaneously incubated.

After washing, the enzyme (streptavidin-peroxydase) is added. After incubation and washing to remove all the unbound enzyme, a substrate solution which is acting on the bound enzyme is added to induce a coloured reaction product. The intensity of this coloured product is directly proportional to the concentration of sCD124 present in the samples.

REAGENTS PROVIDED AND RECONSTITUTION

REAGENTS (Store at 2-8°C)	QUANTITY	COLOUR CODE	RECONSTITUTION
96-wells microtiter plates	2		Ready-to-use
Plastic cover	4		
Standard : 1000 pg/ml	4 Vials Freeze Dried	Yellow	0.51 ml of adequate standard buffer diluent for 1 vial
Standard Buffer Diluent	1 Vial (25 ml)	Black	Dilute in distilled Water
Biotinylated anti-CD124	2 Vials (0.3 ml)	Red	Dilute in biotinylated antibody diluent
Biotinylated Antibody Diluent	1 Vial (13 ml)	Red	Ready-to-use
Streptavidin-HRP	4 Vials (5µl)		0.5 ml of HRP-Diluent before further dilutions
HRP Diluent	1 Vial (23 ml)	Red	Ready-to-use
Washing Buffer	2 Vials (10 ml)	White	Dilute in distilled Water
Chromogen TMB :	1 Vial (21 ml)		Ready-to-use
H2SO4 : Stop Reagent	2 Vials (11 ml)	Black	Ready-to-use

MATERIAL REQUIRED BUT NOT PROVIDED

- * Distilled water.
- * Pipettes : 10 µl, 50 µl, 100 µl, 200 µl and 1000 µl.
- * Vortex mixer and magnetic stirrer.



Cell Sciences, Inc.
480 Neponset Street
Bldg 12A
Canton, MA 02021

Toll Free: 888-769-1246
Phone: 781-828-0610
Fax: 781-828-0542

E-mail: info@cellsciences.com
Web Site: www.cellsciences.com

SAFETY

- * For research use only.
- * The human blood components included in this kit have been tested and found non reactive for HBsAg and anti-HIV. Nevertheless, no known method can offer complete assurance that human blood derivatives will not transmit hepatitis, AIDS or other infections. Therefore, handling of reagents, serum or plasma specimens should be in accordance with local safety procedures , e.g. CDC/NIH Health manual : " Biosafety in Microbiological and Biomedical Laboratories" 1984.
 - * Avoid any skin contact with H₂SO₄ and TMB. In case of contact, wash thoroughly with water.
 - * Do not eat, drink, smoke or apply cosmetics where kit reagents are used.
 - * Do not pipette with mouth.

PROCEDURAL NOTES/LAB. QUALITY CONTROL

1. When not in use, kit components should be stored refrigerated or frozen as indicated on vials or bottles labels. All reagents should be warmed to room temperature before use. Lyophilized standards should be discarded right after resuspension and use.
2. Once the desired number of strips has been removed, immediately reseal the bag to protect the remaining strips from deterioration.
3. Cover or cap all reagents when not in use.
4. Do not mix or interchange reagents between different kit lots.
5. Do not use reagents beyond the expiration date of the kit.
6. Use a clean disposable plastic pipette tip for each reagent, standard, or specimen addition in order to avoid cross-contamination ; for the dispensing of H₂SO₄ and substrate solution, avoid pipettes with metal parts.
7. Use a clean plastic container to prepare the washing solution.
8. Thoroughly mix the reagents and samples before use by agitation or swirling.
9. All residual washing liquid must be drained from the wells by efficient aspiration or by decantation followed by tapping the plate forcefully on absorbent paper. Never insert absorbent paper directly into the wells.
10. The TMB solution is light sensitive. Avoid prolonged exposure to light. Also, avoid contact of the TMB solution with metal to prevent colour development. Warning TMB is toxic avoid direct contact with hands. Dispose off properly.
11. Lyophilised standard vials contain Phenylmethyl Sulfonyl Fluoride (PMSF) which is very toxic. Handle with gloves.
12. If a dark blue colour develops within a few minutes after preparation, this indicates that the TMB solution has been contaminated and must be discarded. Read absorbances within 1 hour after completion of the assay.
13. Serum and plasma samples should be collected in pyrogen/endotoxin-free tubes.
14. Samples should be stored frozen if not analyzed shortly after collection. Avoid multiple freeze-thaw cycles of frozen specimens. Thaw completely prior to analysis.
15. When possible, avoid use of badly hemolyzed or lipemic sera. If large amounts of particles are present, this should be removed prior to assay by centrifugation or filtration.
16. When pipetting reagents, maintain a consistent order of addition from well-to-well. This will ensure equal incubation times for all wells.
17. Respect incubation times described in the assay procedure.
18. Dispense the TMB solution within 15 min. following the washing of the microtiter plate.



REAGENTS PREPARATIONS

Standard buffer diluent

Dilute 5 ml of standard buffer diluent in 50 ml of distilled water or all the content of the vial (25 ml) in 250 ml of distilled water.

Standards

Standards have to be reconstituted with 0.51ml of standard buffer diluent to give a concentration of 1000 pg/ml sCD124.

Dilution of biotinylated anti-CD124

Dilute the biotinylated anti-CD124 with the biotinylated antibody diluent in a clean glass vial according to the number of wells to be used. See the next table for volumes to pipette. Extemporaneous preparations are recommended.

Number of Wells used	Biotinylated Antibody (µl)	Biotinylated Antibody Diluent (µl)
16	40	1060
24	60	1590
32	80	2120
48	120	3180
96	240	6360

Dilution of Streptavidin-HRP

Add 0.5 ml of HRP diluent to a 5 µl vial of Streptavidin-HRP . DO NOT KEEP THIS DILUTION FOR FURTHER EXPERIMENTS. Extemporaneous preparations are recommended. Following the number of wells to be used, further dilutions of Streptavidin-HRP should be made with HRP diluent in a clean glass vial : see hereafter the table for volumes to pipette.

Number of Wells	Streptavidin-HRP(µl)	Strep-HRP Diluent (ml)
16	30	2
24	45	3
32	60	4
48	75	5
96	150	10

Dilution of Washing Buffer

Dilute 2 ml of washing buffer in 400 ml of distilled water or all the content of one vial (10 ml) in 2000 ml of distilled water.



ASSAY METHOD

- a) Before use, mix all reagents thoroughly without making foam.
- b) Determine the number of microwell strips required to test the desired number of samples, plus appropriate number of wells needed for running blanks and standards. Each sample, standard, blank and optional control samples should be assayed in duplicate. Remove sufficient microwell strips from the pouch.
- c) Add 100 µl of appropriate standard diluent to standard wells B1, B2, C1, C2, D1, D2, E1, E2, F1, F2. Reconstitute standard vial with the appropriate volume as described in the chapter reagents preparation. Pipet 200 µl of standard into wells A1 and A2 (see Plate Scheme below). Transfer 100 µl from A1 and A2 to B1 and B2 wells. Mix the contents by repeated aspirations and ejections. Take care not to scratch the inner surface of microwells. Repeat this procedure from the wells B1, B2 to wells C1, C2 and from wells C1, C2 to D1, D2 and so on creating two parallel rows of sCD124 standard dilutions ranging from 1000 to 31.25 pg/ml. Discard 100 µl from the content of the last microwells used (F1, F2).
- d) Add 100 µl of appropriate standard diluent to the blank wells (G1-G2).
- e) Add 100 µl of samples to sample wells.
- f) Preparation of biotinylated anti-CD124 : see reagents preparation.
- g) Add 50 µl of diluted biotinylated anti-CD124 to all wells.
- h) Cover with a plate cover and incubate for 1 hour at room temperature (18°C - 25°C).
- i) Remove the cover plate and wash the plate as follows:
 - 1) aspirate the liquid from each well ;
 - 2) dispense 0.3 ml of washing solution into each well ;
 - 3) aspirate again the content of each well ;
 - 4) Repeat steps 2) and 3) two times.
- j) Prepare HRP solution just before use: see reagents preparation.
- k) Dispense 100 µl of HRP solution into all wells, including the blank wells. Put back the cover plate.
- l) Incubate the microwell strips at room temperature for 30 minutes.
- m) Remove plate cover and empty wells. Wash microwell strips according to point (i). Proceed immediately to the next step.
- n) Pipette 100 µl of ready-to-use TMB substrate solution into all wells, including the blank wells and incubate in the dark for 12-15 minutes at room temperature. Avoid direct exposure to light by wrapping the plate in aluminium foil.
- o) Incubation time of the substrate solution is usually determined by the ELISA reader performances: many ELISA readers record absorbance only up to 2.0 O.D. The O.D. values of the plate should be monitored and the substrate reaction stopped before positive wells are no longer properly readable (maximum 20 minutes).
- p) The enzyme-substrate reaction is stopped by quickly pipetting 100 µl of 1.8 N sulfuric acid into each well, including the blank wells, to completely and uniformly inactivate the enzyme. Results must be read immediately after the addition of sulfuric acid, or within one hour, if the microwell strips are stored at 2-8°C in the dark.
- q) Read absorbance of each well on a spectrophotometer using 450 nm as the primary wavelength and optionally 620 nm (610 nm to 650 nm is acceptable) as the reference wavelength.




SUGGESTED PLATE SCHEME

	Standard Concentrations pg/mL		Sample wells									
	1	2	3	4	5	6	7	8	9	10	11	12
A	100 0	100 0										
B	500	500										
C	250	250										
D	125	125										
E	62.5	62.5										
F	31.2 5	31.2 5										
G	Blan k	Blan k										
H												

DATA ANALYSIS

Generate a linear standard curve by plotting the average absorbance on the vertical axis versus the corresponding sCD124 standard concentration on the horizontal axis. The amount of sCD124 in each sample is determined by extrapolating OD values to sCD124 concentrations using the standard curve.

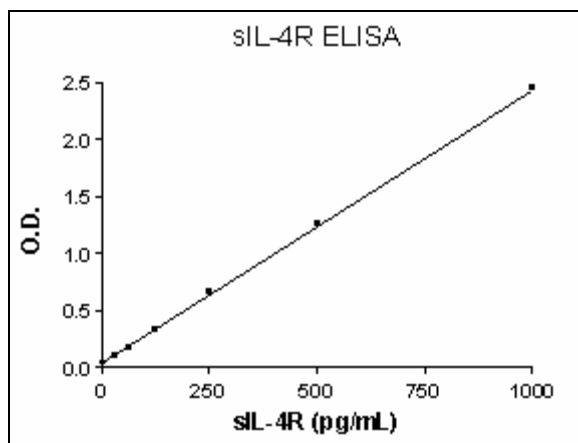

Typical sCD124 standard curve ranging from 31.25 to 1000 pg/mL



LIMITATIONS OF THE PROCEDURE

Do not extrapolate the standard curve beyond the 1000 pg/ml standard curve point. The dose-response is non-linear in this region and good accuracy is difficult to obtain. Concentrated samples (> 1000 pg/ml) have to be diluted with standard diluent or with your own sample buffer. During analysis, multiply results by the appropriate dilution factor.

The influence of various drugs, aberrant sera (hemolyzed, hyperlipidemic, jaundiced...) and the use of plasma instead of serum specimens has not been thoroughly investigated. The rate of degradation of native sCD124 in various matrices has not been investigated. The immunoassay literature contains frequent references to aberrant signals observed with some human sera and attributed to heterophilic (or anti-mouse) antibodies. Though we have not found such specimens to date, the possibility of this occurrence can not be excluded.

PERFORMANCES AND CHARACTERISTICS

Sensitivity

The minimum detectable dose of sCD124 is less than 11.5 pg/ml. This has been determined by adding 3 standard deviations to the mean optical density obtained when the zero standard was assayed 32 times.

Precision

Intra-Assay					Inter-Assay				
Sample	n	Mean (pg/mL)	SD	CV%	Sample	n	Mean (pg/mL)	SD	CV%
A	8	563.7	17.8	3.1	A	8	530.5	35.1	6.6
B	8	140.6	9.6	6.8	B	8	106.5	7.58	7.1



ASSAY PROCEDURE SUMMARY Total procedure length : 1h45

Add 100µl of sample or diluted standard

|

Add 50µl of diluted biotinylated
detection antibody to all wells

|

Incubate 1 hour at room temperature

|

Add 100µl of HRP-Streptavidin conjugate

|

Incubate 30min at room temperature

|

Add 100 µl of ready-to-use TMB
Protect from light. Let the color develop for 10-15 mn.

|

Add 100µl H₂SO₄

|

Read Absorbance at 450nm



REFERENCES

Fanslow, WC. et al. (1990) : A soluble form of the interleukin 4 receptor in biological fluids. Cytokines 2 : 398-401.

Keegan, AD. Et al. (1991) : The IL-4 receptor : biochemical characterization of IL-4-binding molecules in a T cell line expressing large numbers of receptors. J. Immunol. 146 : 2272-2279.

Fanslow, WC et al. (1991) : Regulation of alloreactivity in vivo by IL-4 and the soluble IL-4 receptor. J. Immunol. 147 : 535-540.

Fanslow, WC et al. (1993) : Identification of a distinct low-affinity receptor for human interleukin-4 on pre-B cells. Blood. 81 : 2998-3005.

Breit, S. et al. (1996) : A strict requirement of interleukin-4 for interleukin-4 induction in antigen-stimulated human memory T cells. Eur. J. Immunol. 26 : 1860-1865.

Sang, DK. (1999) : Increased levels of soluble interleukin-4 receptor in the sera of patients with visceral leishmaniasis. J. Infect. Dis. 179 : 743-746.

Jung, T. (1999) : Enhancement of human IL-4 activity by soluble IL-4 receptors in vitro. Eur. J. Immunol. 29 : 864-871.

